

MIE 402 Lab 2

Single Pendulum | Fall 2026 | Prof. Ehsan Roohi

Your group will record a single pendulum with a high-speed camera, use Pendulum Motion Analyzer to track the bob, convert tracked position to angle, and compare the measured response with a theoretical model. The three trials test the limit of the small-angle approximation.

Learning objectives

- Identify the physical meaning of pivot, bob, length, angle, and angular velocity.
- Use high-speed video and image tracking for dynamic measurement.
- Convert bob coordinates into angular displacement about a fixed pivot.
- Compare measured period and frequency with theory.
- Explain why finite initial angle changes pendulum period.

Equipment and required trials

- Single-pendulum assembly, black background, protractor, calibration marks, and scale.
- High-speed camera, tripod, lens, ring light, power supply, and SD card.
- Computer with MATLAB and Pendulum Motion Analyzer.
- Trial A: less than 10 degrees; Trial B: 10 to 90 degrees; Trial C: 90 to 180 degrees.
- Every trial begins from rest. Record initial angle, length to bob center, relevant masses, frame rate, calibration distance, and file name.



Figure 1. Single-pendulum apparatus. Track the center of the bob and use the pivot rotation center as reference.

Camera setup and safety

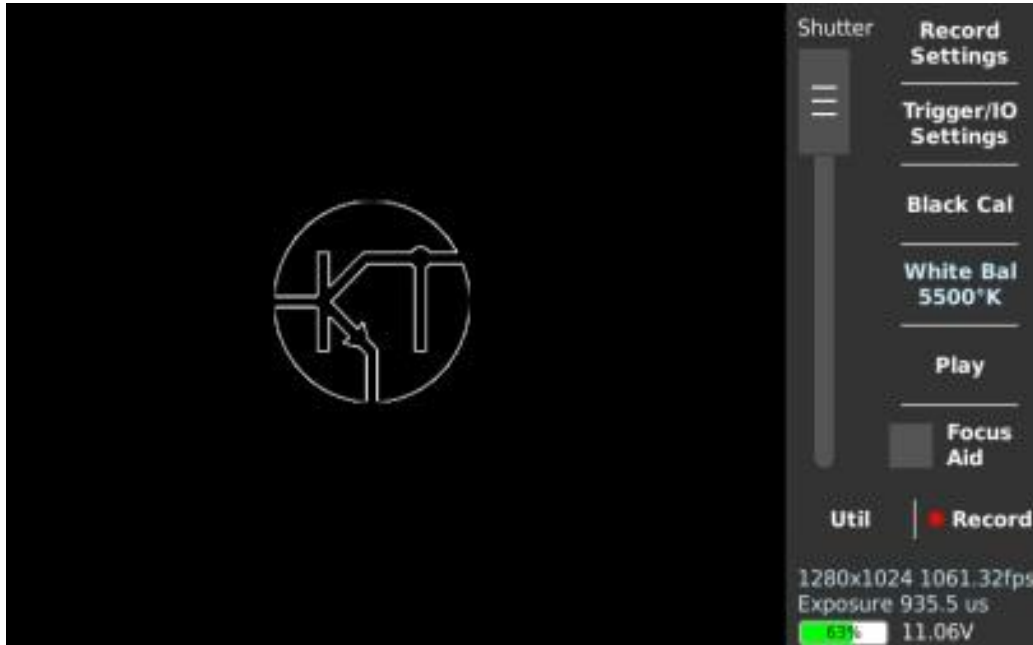


Figure 2. Camera main screen. Record, playback, shutter, and Record Settings are available here.

1. Mount camera and lens securely; mount ring light as shown by the GTA.
2. Set pendulum against black background. Confirm pivot and bob have high contrast.
3. Insert SD card gently, connect power, and confirm tripod stability.
4. Frame the full motion and perform a low-amplitude test swing.
5. Adjust focus and brightness until the bob edge is sharp without reflection at its center.

Walk carefully around camera, tripod, and cables. Do not move the tripod while recording.

Camera settings before each trial

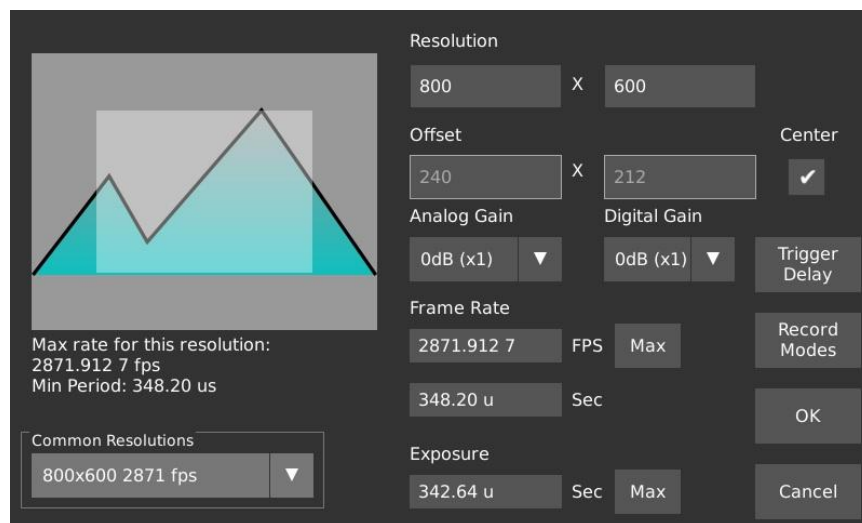


Figure 3. Record Settings. Use the GTA-assigned resolution and frame rate. The normal target is 1280 x 1024 at 1000 frames/s when available.

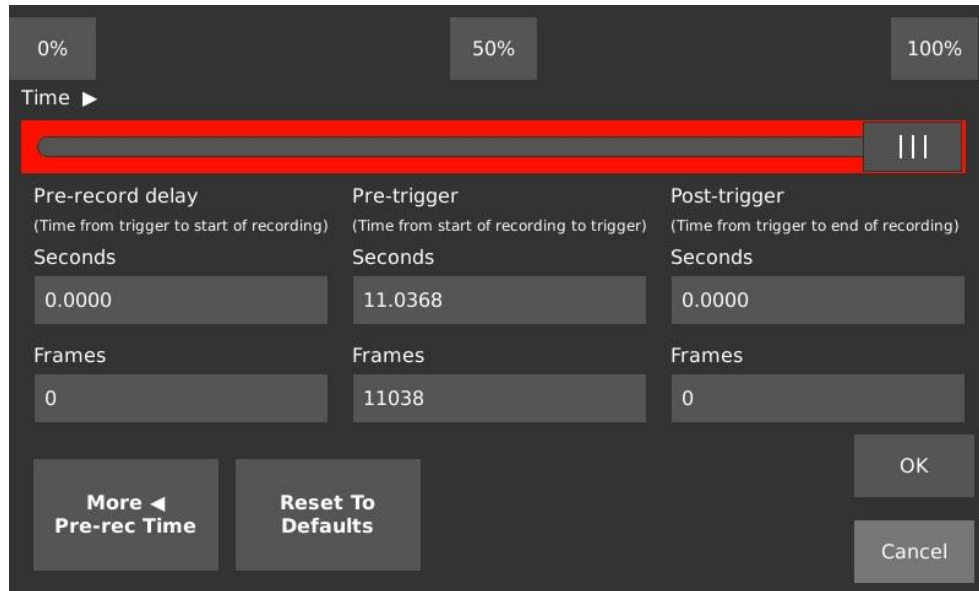


Figure 4. Trigger delay. Choose the GTA-specified setting so release and several full cycles appear in the saved video.

- Record actual resolution and frame rate in notes.
- Use trigger delay that captures release and several cycles.
- Use test video to check focus, contrast, and full motion range.
- Name every record with angle and trial, for example single_35deg_trial01.

Record and save

6. Press Record.
7. Hold pendulum at measured initial angle, keep it stationary, and release without pushing.
8. After about 15 seconds, press Stop.
9. Use Playback to confirm release and early cycles are visible.
10. Mark useful segment. Select H.264, turn Text Overlay off, save to SD card, and wait for completion.
11. Safely remove SD card and copy videos to group storage.

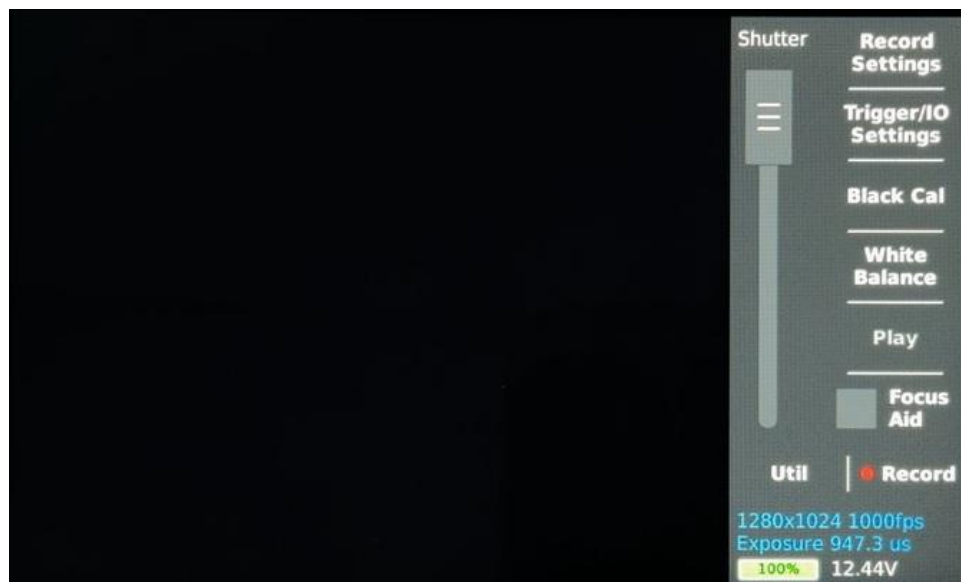


Figure 5. Save settings. Confirm filename, H.264 format, Text Overlay off, and SD-card location.

Track in Pendulum Motion Analyzer

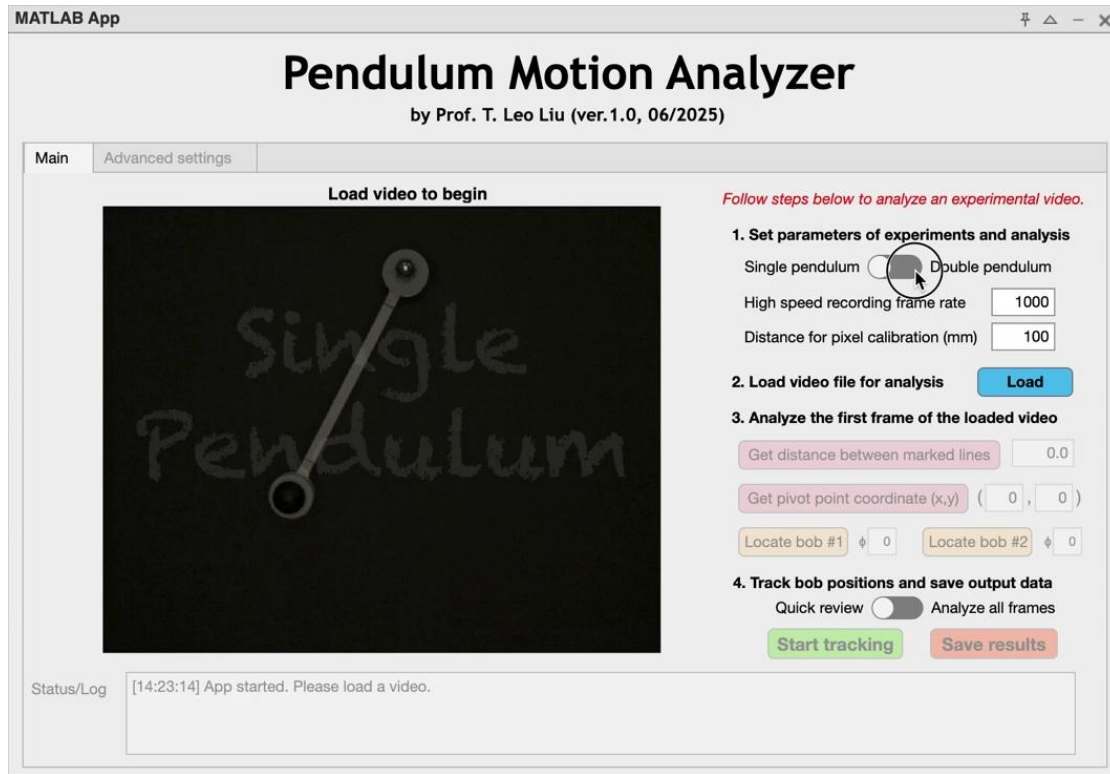


Figure 6. Select Single pendulum and enter the actual camera frame rate.

12. Enter camera frame rate and 100 mm calibration distance between black arm marks unless the GTA says otherwise.
13. Load video, draw calibration line in first usable frame, then locate pivot and bob center.
14. Run Quick Review. If the marker misses more than about 10 percent of frames, improve focus or repeat point selection.
15. Run Analyze All Frames after review succeeds. Save Results as MAT file.

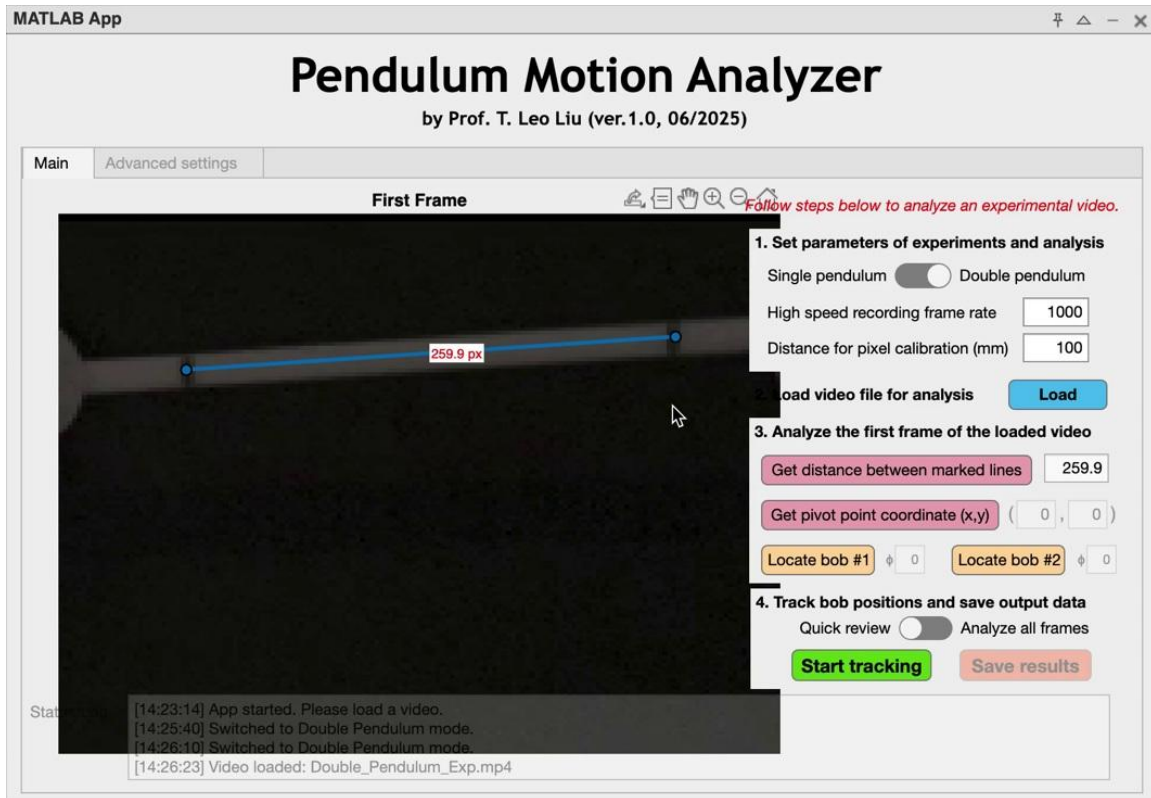


Figure 7. The known distance between black marks converts pixels to millimeters.

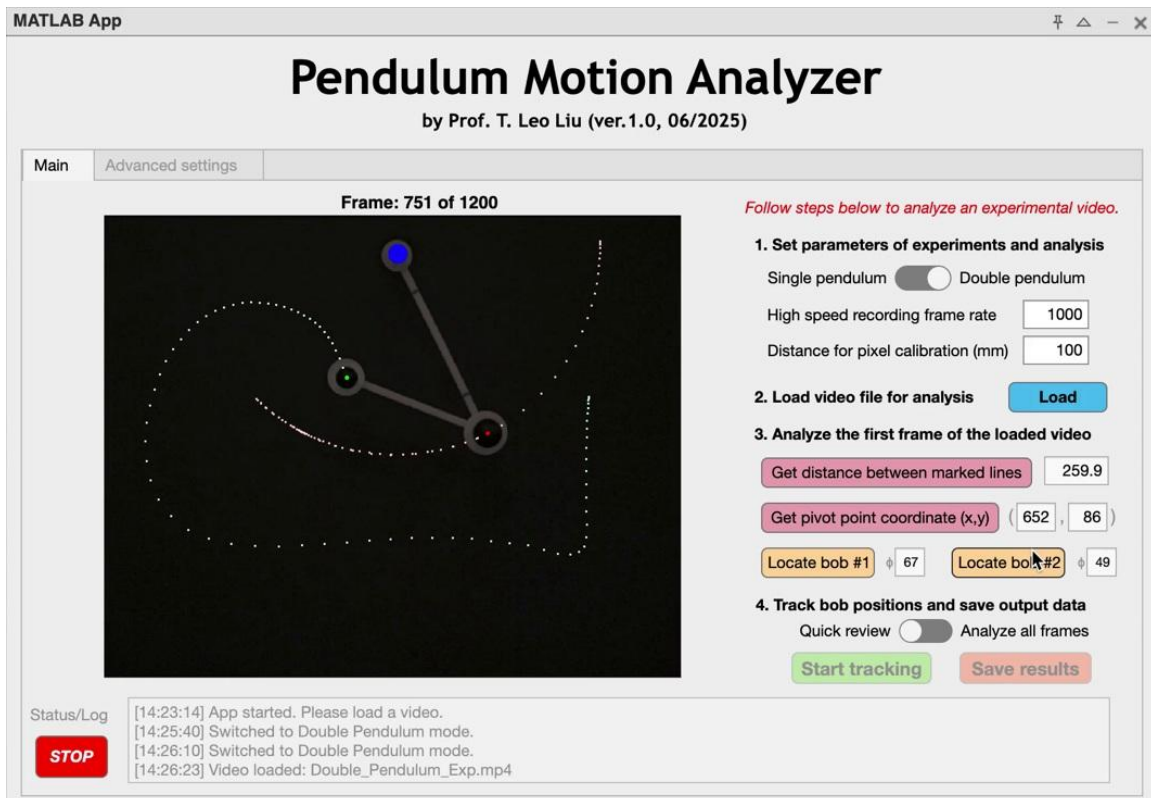


Figure 8. The overlay marker should follow the bob center through the usable record.

MAT output and analysis

The saved ExpData matrix contains time and tracked coordinates. For a single pendulum use time and bob-1 x and y; confirm exact column order in the app. With pivot x_p , y_p and bob x_b , y_b , calculate $\theta = \text{atan2}(x_b - x_p, y_b - y_p)$. Use radians in equations and label degree plots clearly.

- For each trial plot $\theta(t)$ and angular velocity. State derivative method and any smoothing.
- Plot spectra for θ and angular velocity. Report dominant frequency and FFT-bin spacing.
- Calculate small-angle frequency $\sqrt{g/L}/(2\pi)$. For larger angles compare with nonlinear model or explain small-angle limitation.
- Create a table of length, mass, initial angle, theoretical frequency, time-domain frequency, spectral frequency, and meaningful differences.
- Discuss camera calibration, tracking error, air resistance, pivot friction, finite angle, and release error.

Before leaving

- One saved video and verified MAT file exist for each angle range.
- Notes include frame rate, calibration distance, length, masses, initial angle, and file names.
- Group checks one $\theta(t)$ plot and one spectrum with GTA.
- Back up raw videos, MAT files, notes, and exported figures.